


## Miniature Circular ESF Components and Stater Kit

PART NO. DESCRIPTION

4.70007 Spherical Nut, 4mm

4.70008 Spherical Washer, 4mm

By placing a pair on each side of a ring, rod angle can be adjusted to fine-tune reduction and assembly. Six pairs will allow a ring or ring block with three rods to be adjusted. **Must be used as mating pairs** requiring two 7mm wrenches. Stainless steel.

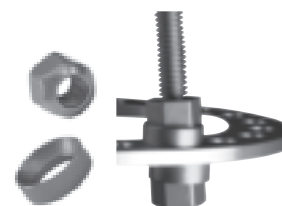

300.003 Hex Bolt, 4mm x 8mm

Used as a general-purpose fastener with the mini CESF system. This bolt is also a component of the mini SK® ESF clamp assembly and may be useful for drop wire technique off rings or off one-hole posts. In addition, this bolt can be used to mount mini angular hinges to rings but may require flat washers acting as spacers. Stainless steel.



### Mini Circular ESF Starter Kit

Part No. 4.00001

Qty.	Part No.	Description	Qty.	Part No.	Description
4	4.13501	Full Ring, 35mm x 2.4mm	18	4.30001	Cannulated/Slotted Wire
3	4.14501	Full Ring, 45mm x 2.4mm			Fixation Bolt
1	4.14504	Stretch Ring, 45mm x 2.4mm	6	4.30005	Slotted Wire Fixation Washer
3	4.26050	Threaded Rod, 4mm x 50mm	6	300.003	Hex Bolt, 4mm x 8mm
3	4.26075	Threaded Rod, 4mm x 75mm	50	4.30006	Hex Nut, 4mm
3	4.26100	Threaded Rod, 4mm x 100mm	12	4.30008	Flat Washer, 4mm
1	4.50005	7mm Double Wrench	2	4.50003	7mm Combination Wrench
8	4.70001	Nylon Zero Tolerance Nut, 4mm	1	4.50004	7mm Angled Socket Wrench
2	4.70005	Mini One-hole Post	8	4.40002	Stopper Wire, 1.0mm x 125mm